

PRODUCT INFORMATION

HUMAN iPSC-DERIVED DOPAMINERGIC NEURON

Through our Human iPSC-derived dopaminergic neurons you have the potential to gather meaningful insights into new treatments, diagnostics, and disease prevention.

Human iPSC-Derived Dopaminergic Neurons

Catalog Number	40HU-002-1M 40HU-002-2M	Cell Number	1.0 million cells/vial (Cryopreserved) 2.0 million cells/vial (Cryopreserved)
Species	<i>Homo sapiens</i>	Storage Temperature	Liquid nitrogen

Product Description

Midbrain dopamine (DA) neurons are critical for directing fundamental brain functions such as voluntary movement, reward processing, and working memory. The substantia nigra (SN) and the ventral tegmental area (VTA) have the highest populations of DA neurons in the midbrain. The degeneration of DA neurons within the pars compacta region of the SN is a pathological hallmark of Parkinson's disease (PD) and Lewy body dementia (LBD)^[1]. For many years, powerful experimental model organisms like the mouse, fruit fly, and baker's yeast have been used to study neurodegenerative diseases, providing insights into disease mechanisms like pathological aggregation of key proteins, the nature and processes of neuronal damage, the role of genetic determinants, and the contribution of neuroinflammation in fueling neuronal loss^{[2][3]}. However, it appears that the use of these models has only partially elucidated some elements of the illnesses, impeding a meaningful translation into new treatments, diagnostics, and prevention.

iXCells Biotechnologies is proud to provide fully differentiated and functional human iPSC-derived DA neurons that display typical neuronal morphology and express all key markers of DA neurons, e.g., TH, FoxA2 (Figure 1) when cultured in the Human Dopaminergic Neuron Maturation Medium (Cat# MD-0105-100ML). In addition, our iPSC-derived DA neurons can also be co-cultured with glial cells or other cell types for drug screening platforms.

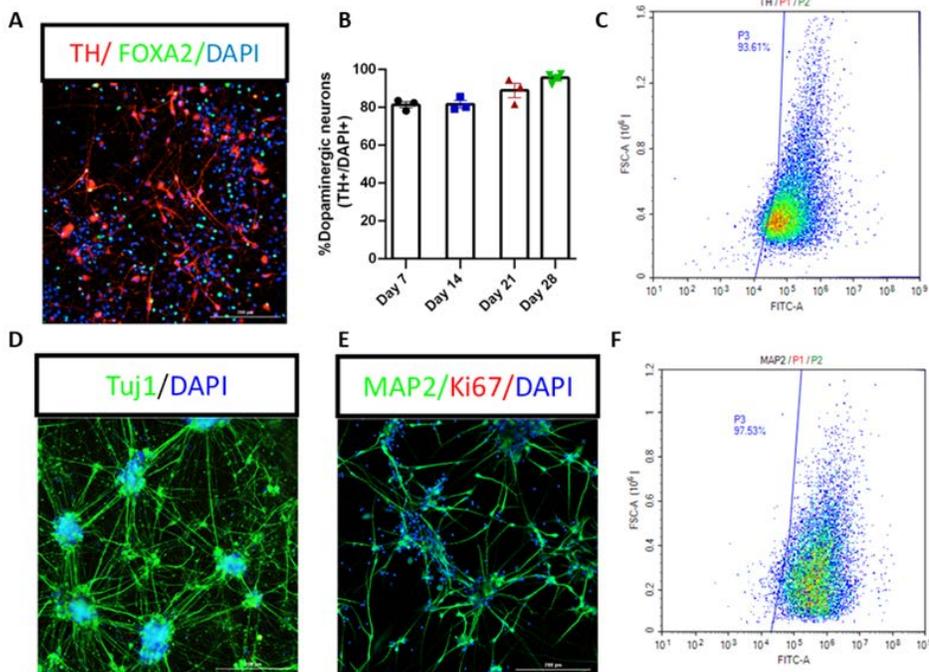


Figure 1. Human iPSCs derived dopaminergic neurons show expression of characteristic biological markers. (A) Immunostaining shows the expression of midbrain dopaminergic neuron markers FoxA2 and TH, 28 days post-thawing. (B) Quantifications for percentage of dopaminergic neurons, positive for the TH marker over multiple time-points. Results are expressed as means \pm SEM. (C) Flow cytometry measurements demonstrate a highly specific population of fully differentiated midbrain dopaminergic neuron (TH). (D)(E) Immunostaining shows the expression of neuron marker Tuj1 (D) 21 days post-thawing and MAP2 (E) 7 days post-thawing and the absence of the proliferative progenitor marker Ki67. (F) Flow cytometry measurements demonstrate a highly pure population of fully differentiated neurons (MAP2). Nuclei were counterstained with DAPI. Scale bars, 200 μ m.

Product Details

Tissue Origin	Human iPSC-derived dopaminergic neurons (Normal)
Package Size	1.0 million cells/vial; 2.0 million cells/vial;
Shipped	Cryopreserved
Storage	Liquid Nitrogen
Media	Human Dopaminergic Neuron Maturation Medium (Cat # MD-0105-100ML)

Protocols - Monoculture of hiPSC-Derived Dopaminergic Neurons

The following protocol is based on 48-well plate format

1. Prepare the coating vessel before thawing the cells. Coated plate preparation: coat the plates with Poly-L-ornithine solution (10 µg/ml) in water overnight in the incubator, wash twice with water the next day and coat the plates with laminin (5 µg/ml) in dPBS overnight. Do not let dry at any stage.

Note: Upon receipt of the frozen cells, it is recommended to thaw the cells and initiate the culture immediately in order to retain the highest cell viability.

2. To thaw the cells, put the vial in 37°C water bath with gentle agitation for ~1-2 minutes. Keep the cap out of water to minimize the risk of contamination.

3. Pipette the cells into a 15 mL conical tube with 5 mL **Human Dopaminergic Neuron Maturation Medium (Cat# MD-0105-100ML)**.

4. Centrifuge at 300g for 5 minutes at room temperature.

5. Remove the supernatant and re-suspend the cells in **Human Dopaminergic Neuron Maturation Medium + 10 µM Y27632**. Y27632 should be removed after 24 hours by replenishing fresh medium carefully.

6. Seed the cells on precoated plate at the desired density. Incubate in 37°C CO₂ incubator overnight.

Note: We recommend seeding 90-200K cells per cm² depending on the application. Cell debris may be observed after cell recovery. Refer to the CoA of each lot to determine the seeding density for your experiment.

7. Perform half medium change every 2-3 days. Most of the cells should express high levels of Tuj1 and MAP2, 2 days after thaw, and express high levels of TH and FOXA2 7-10 days after thaw.

Note: DA neurons tend to aggregate and detach from the plates. Change 50% of the medium with extra care to avoid cell loss.

REFERENCE

1. Bloem, B. R., Okun, M. S. & Klein, C. (2021) Parkinson's disease. *Lancet* 397, 2284–2303.
2. Schulz-Schaeffer, W.J. (2015). Is Cell Death Primary or Secondary in the Pathophysiology of Idiopathic PD? *Biomolecules* 5, 1467-1479.
3. Poewe, W., Seppi, K., Tanner, C. et al. (2017). Parkinson disease. *Nat Rev Dis Primers* 3, 17013

Find out more at [ixcellsbio.com/product](https://www.ixcellsbio.com/product)

This product is intended for laboratory research purposes only. It is not intended for use in humans. While iXCells Biotechnologies uses reasonable efforts to include accurate and up-to-date information on this product sheet, we make no warranties or representations as to its accuracy. Citations from scientific literature and patents are provided for informational purposes only. iXCells Biotechnologies does not warrant that such information has been confirmed to be accurate. This product is sent with the condition that you are responsible for its safe storage, handling, and use. iXCells Biotechnologies is not liable for any damages or injuries arising from receipt and/or use of this product. While reasonable effort is made to insure authenticity and reliability of strains on deposit, iXCells Biotechnologies is not liable for damages arising from the misidentification or misrepresentation of cultures.



Address: 7270 Trade Street, Suite 102,
San Diego, CA 921216

Orders: orders@ixcellsbiotech.com

Technical Support : supports@ixcellsbiotech.com

Phone: +1 (858)412-5988