

## Product Information


### Human Sertoli Cells (HSerC)

Catalog Number	10HU-149	Cell Number	1.5 million cells/vial
Species	<i>Homo sapiens</i>	Storage Temperature	Liquid nitrogen

### Product Description

Human Sertoli Cells (HSerC) are essential for testicular development, spermatogenesis, and formation of the blood-testis barrier [1, 2]. HSerC limit the passage of substances such as hormones and nutrients to the adluminal compartment of the seminiferous tubules [1]. In addition to forming the blood-testis barrier, HSerC also provide the main structural support for the seminiferous tubules and protect the germ cells from the immune system [1]. Aberrant HSerC proliferation can contribute to the development of male reproductive disorders such as testicular germ-cell cancer, cryptorchidism, hypospadias, and low sperm count [2]. HSerC proliferation is in part controlled by follicle-stimulating hormone (FSH) and thyroid hormone (TH), where FSH drives proliferation and TH promotes a more quiescent state [3]. Cultured HSerC are a useful in vitro model to better understand testicular dysgenesis syndrome and to develop treatments for male reproductive disorders.

iXCells Biotechnologies provides high quality HSerC, which are isolated from human testis and cryopreserved at P1, with >1.5 million cells in each vial. HSerC express GATA-4 and Sox-9 and are negative for HIV-1, HBV, HCV, mycoplasma, bacteria, yeast, and fungi. They can further expand for no more than 3 passages in Sertoli Cell Growth Medium (Cat# MD-0091) under the condition suggested by iXCells Biotechnologies. Further expansion may decrease the purity.

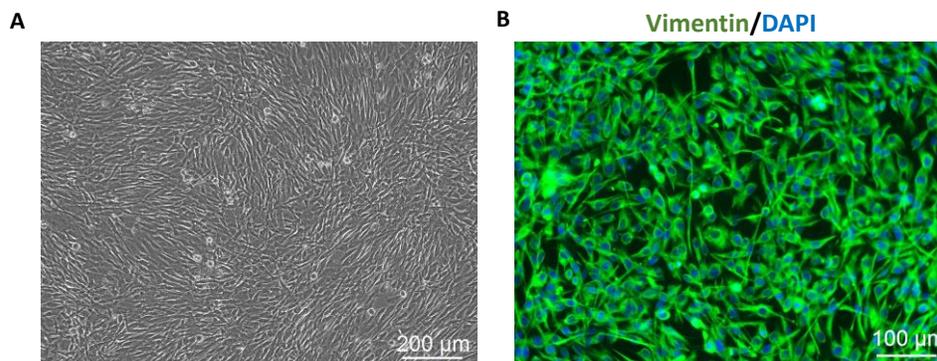

**Figure 1.** Human Sertoli Cells (HSerC). (A) Phase contrast image of HSerC. (B) Immunofluorescence staining with antibody against Vimentin.

## Product Details

<b>Tissue</b>	Human testis
<b>Package Size</b>	1.5 million cells/vial
<b>Passage Number</b>	P1
<b>Shipped</b>	Cryopreserved
<b>Storage</b>	Liquid nitrogen
<b>Growth Properties</b>	Adherent
<b>Media</b>	Sertoli Cell Growth Medium (Cat# MD-0091)

## Protocols

### Thawing of Frozen Cells

1. Upon receipt of the frozen HSerC, it is recommended to thaw the cells and initiate the culture immediately in order to retain the highest cell viability.
2. To thaw the cells, put the vial in 37°C water bath with gentle agitation until the contents completely thaw. Keep the cap out of water to minimize the risk of contamination.
3. Carefully transfer the cells into a 15 mL conical tube with ~5 mL fresh **Sertoli Cell Growth Medium (Cat# MD-0091)**. Gently re-suspend and dispense the contents.
4. Transfer the cells into the desired culture vessels. A seeding density of 5,000 cells/cm<sup>2</sup> is recommended
5. Place the culture vessels to the 37°C incubator at 5% CO<sub>2</sub> for continuous culture.
6. For the best result, do not disturb the culture for at least 16 hours after the culture has been initiated. Refresh culture medium the next day to remove unattached cells, then every other day thereafter.

**Safety Precaution:** *it is highly recommended that protective gloves and clothing should be used when handling frozen vials.*

### Standard Culture Procedure

1. HSerC should be maintained in the **Sertoli Cell Growth Medium (Cat# MD-0091)**.
2. Change the medium every other day until cells reach ~90% confluence.
3. Remove the medium, and wash once with sterile PBS (5 mL/T75 flask).
4. Add ~3 mL of 0.25% Trypsin-EDTA to the flask and incubate for 3-5 minutes at 37°C. Neutralize the enzyme by adding 6-9 mL culture medium.
5. Centrifuge at 1,000 rpm (~220 g) for 5 minutes, and re-suspend the cells in desired volume of culture medium.
6. Transfer the cells to the desired culture vessels with the recommended cell density.

## References

- [1] Chui K, Trivedi A, Cheng C, Cherbavaz, Dazin P, Huynh A, Mitchell J, Rabinovich G, Noble-Haeusslein L, John C. (2011) "Characterization and functionality of proliferative human Sertoli cells." *Cell Transplant*. 20(5): 619-635.
- [2] Sharpe R, McKinnell C, Kivlin C, Fisher J. (2003) "Proliferation and functional maturation of Sertoli cells, and their relevance to disorders of testis function in adulthood." *Reproduction*. 125: 769-784.
- [3] Tarulli G, Stanton P, Meachem S. (2012) "Is the adult Sertoli cell terminally differentiated" *Biol Reprod*. 87(1): 1-11.

## Disclaimers

This product is intended for laboratory research purposes only. It is not intended for use in humans. While iXCells Biotechnologies uses reasonable efforts to include accurate and up-to-date information on this product sheet, we make no warranties or representations as to its accuracy. Citations from scientific literature and patents are provided for informational purposes only. iXCells Biotechnologies does not warrant that such information has been confirmed to be accurate.

This product is sent with the condition that you are responsible for its safe storage, handling, and use. iXCells Biotechnologies is not liable for any damages or injuries arising from receipt and/or use of this product. While reasonable effort is made to insure authenticity and reliability of strains on deposit, iXCells Biotechnologies is not liable for damages arising from the misidentification or misrepresentation of cultures.

© iXCells Biotechnologies 2015. All rights reserved.