

Product Information

Mouse Epidermal Keratinocytes – Adult (MEK-a)

Catalog Number	10MU-035	Cell Number	0.5 million cells/vial
Species	<i>Mus Musculus</i>	Storage Temperature	Liquid nitrogen

Product Description

Epidermal keratinocyte is the predominant cell type in the outermost layer of the skin - the epidermis, which serves as a critical barrier to separate and protect the inside of human body from outside environment and damage from pathogens, heat, UV radiation and water loss. Epidermal keratinocytes originate in the stratum basale, and they undergo gradual differentiation and migrate towards the surface of the epidermis until they reach the stratum corneum, where they form a tight layer of nucleus-free and highly keratinized squamous cells. This layer forms an effective barrier to prevent water loss and the entry of infectious agents. Keratinocytes are also known to produce various growth factors, cytokines, antimicrobial peptides, and complement factors. Therefore, keratinocytes are important for wound healing, inflammation, infection, skin microbiome and immune response.

iXCells Biotechnologies provides high quality primary Mouse Epidermal Keratinocytes-adult (MEK-a), which are isolated from adult mouse tail and cryopreserved at P0, with ≥ 0.5 million cells in each vial. MEK-a are characterized by phalloidin staining and are negative for HIV-1, HBV, HCV, mycoplasma, bacteria, yeast and fungi. **MEK-a are not recommended for expanding or long term cultures since the cells do not proliferate in culture.**

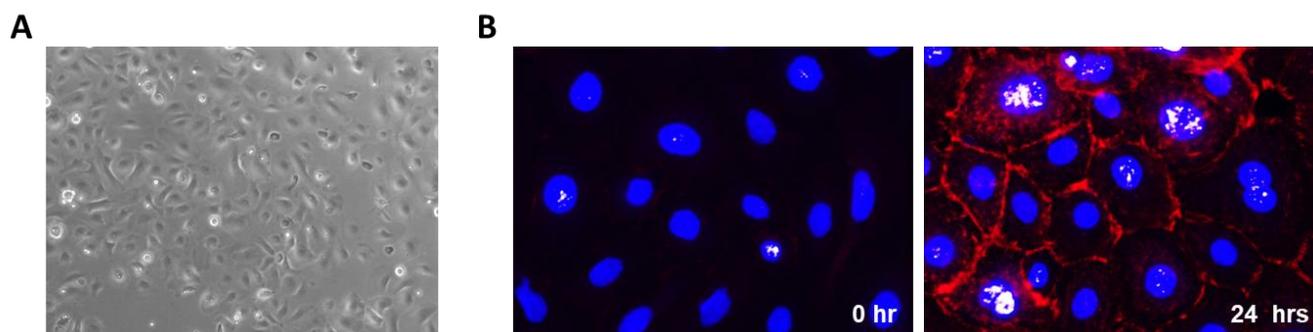


Figure 1. (A). Phase contrast image of primary Mouse Epidermal Keratinocytes - Adult (MEK-a). (B). MEK-a were differentiated in the presence of 0.2mM CaCl_2 , and cells were stained with phalloidin (red) to visualize the formation of actin fiber-rich filopodial projections between adjacent cells during differentiation. Nuclei were counterstained with DAPI.

Product Details

Tissue	Epidermis from adult mouse tail
Package Size	0.5 x10 ⁶ cells/vial (Enough to seed on two wells of 6-well culture plate or four wells of 12-well culture plate, and cells will reach 80% confluency in 4-5 days for assay).
Passage Number	P0
Shipped	Frozen
Storage	Liquid nitrogen
Growth Properties	Adherent
Media	Keratinocyte Growth Medium (Cat# MD-0047)

Protocols

Prepare culture vessels with Coating Matrix Kit (Invitrogen, Catalo# R011K)

1. Dilute 100µl Coating Matrix (50-9700) in 10ml Dilution Medium (50-9701) at 1:100, and mix well;
2. Add 5ml mixed coating matrix solution per each T75 flask;
3. Rock back and forth to ensure uniform distribution of the coating matrix over the surface of the flask;
4. Cap the flasks and incubate for 30 minutes at room temperature;
5. Remove excess Coating Matrix/Dilution Medium from each flask. The flasks may be used immediately, or may be stored at 2° to 8° C for short periods.

Note: Coating the culture vessels is required.

Thawing of Frozen cells and standard culture procedure

1. Upon receipt of the frozen cells, it is recommended to thaw the cells and initiate the culture immediately in order to retain the highest cell viability.
2. To thaw the cells, put the vial in 37°C water bath with gentle agitation for 1-2 minutes. Keep the cap out of water to minimize the risk of contamination.
3. Pipette the cells into a 15 mL conical tube with 5 mL fresh Keratinocyte Growth Medium (Cat# MD-0047).
4. Centrifuge at 1,000 rpm (~220 g) for 5 minutes under room temperature.
5. Remove the supernatant and resuspend the cells in fresh culture medium.
6. Culture the cell in T75 flask precoated with Attachment Factor (Thermo Fisher, Cat#S006100) in a 37°C CO₂ incubator. For best results, do not disturb the culture for at least 24 hours after the culture has been initiated.
7. Change the culture medium every other day until the culture is approximately 50% confluent. Then change the medium daily until the cells reach the desired confluency prior to the experiment.

Note: Upon thawing of the cells, please plate cells as soon as possible because prolong incubation in suspension will lead to terminal differentiation of keratinocytes.

Safety Precaution: it is highly recommended that protective gloves and clothing should be used when handling frozen vials.

References

- [1] Raja, Sivamani K, Garcia MS, Isseroff RR. *Front Biosci.* 2007;12:2849-68. Wound re-epithelialization: modulating keratinocyte migration in wound healing.
- [2] Proksch E1, Brandner JM, Jensen JM. *Exp Dermatol.* 2008 Dec;17(12):1063-72. The skin: an indispensable barrier.
- [3] Pasparakis M, Haase I, Nestle FO. *Nat Rev Immunol.* 2014 May;14(5):289-301. Mechanisms regulating skin immunity and inflammation.

Disclaimers

This product is intended for laboratory research purposes only. It is not intended for use in humans. While iXCells Biotechnologies uses reasonable efforts to include accurate and up-to-date information on this product sheet, we make no warranties or representations as to its accuracy. Citations from scientific literature and patents are provided for informational purposes only. iXCells Biotechnologies does not warrant that such information has been confirmed to be accurate.

This product is sent with the condition that you are responsible for its safe storage, handling, and use. iXCells Biotechnologies is not liable for any damages or injuries arising from receipt and/or use of this product. While reasonable effort is made to insure authenticity and reliability of strains on deposit, iXCells Biotechnologies is not liable for damages arising from the misidentification or misrepresentation of cultures.
© iXCells Biotechnologies 2015. All rights reserved.